




Thermo Scientific Heraeus® Cryofuge® 5500i

**Flexible refrigerated
centrifuges for blood banking
and clinical applications**

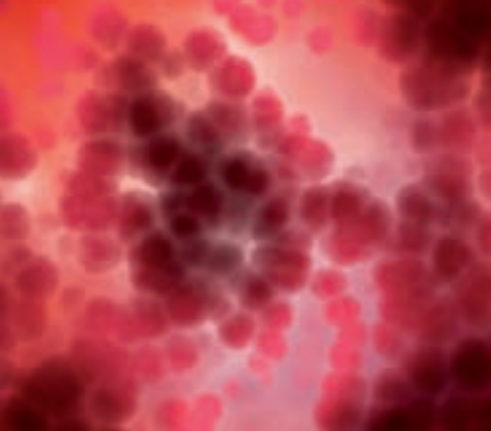

supplyLAB

www.supplylab.pt
geral@supplylab.pt

Cacém Park - Edifício 9
Estrada de Paço de Arcos nº88
2739-512 Agualva Cacém
T +(351) 21 4278700
F +(351) 21 4278709

Thermo Scientific Heraeus Cryofuge 5500i – Refrigerated Blood Processing Centrifuges

The Thermo Scientific Heraeus Cryofuge 5500i offers optimal capacity for small to medium size blood processing centers. The Heraeus Cryofuge 5500i can spin up to 8 blood bags, in a wide range of bag configurations – including 500 mL quints with in-line filters. In addition, achieve greater throughput from higher g-forces and capacities for all popular blood collection types and clinical tubes.

- Optimal 4 or 8 blood bag capacity for component processing
- Thermo Scientific Diagnostik™ rotor spins clinical analyzer racks pre-loaded with tubes, minimizing sample handling steps
- Superior 32 standard microplate capacity
- Thermo Scientific HIGHplate® rotor

- spins microplates, or filter plates up to 85 mm in height, at 5,890 x g
- Easy-to-clean soft-touch keypad
- Reduced maintenance with brushless induction drive

In keeping with the latest trends in clinical testing, the Heraeus Cryofuge 5500i maximizes microplate processing with the added benefit of biocontainment.

A wide range of rotors – including windshielded swinging bucket, clinical analyzer rack and windshielded microplate varieties – offer maximum flexibility for blood component processing and clinical testing needs.



supplyLAB
www.supplylab.pt

Its powerful induction drive and its swinging bucket rotor – which accommodates DoubleSpin™ buckets – combine to provide higher g-forces and capacities for every type of sample tube. The result: shorter processing times for more samples.

Safe operation for you and your sample

The Heraeus Cryofuge 5500i meets stringent international safety standards. It offers the widest

selection of independently tested biocontainment rotors that comply with rigorous international guidelines for aerosol containment of bacteria and viruses.

SMARTspin Imbalance Detection System

provides optimum run stability for all rotors at all speeds, and allows sample tubes to be balanced by eye alone. The system will automatically decelerate the centrifuge if an excessive imbalance is detected.

Automatic Rotor Recognition offers convenient operation as well as reliable protection against incorrect parameter settings, and ensures reproducible run conditions.

Easy instrument set-up and operation includes:

- 9 program memory for standard runs
- 9 acceleration and deceleration profiles
- rpm to RCF conversion
- 'Pretemp' function for fast pre-cooling
- Quick Run Mode

Component processing made easy

The Heraeus Cryofuge 5500i offers convenient operation at your fingertips. Unlike other models, the Thermo Scientific EASYcontrol™ II Soft-Touch Keypad is located at the front of the instrument, eliminating awkward reaching. It incorporates large LED display panels that can be cleaned easily and read even at a distance.

Utilize your existing adapter inventory

The Heraeus Cryofuge 5500i accepts the most popular liners, adapters and accessories from the Heraeus Cryofuge 6000i and/or Heraeus Multifuge 3+:

- All accessories available for the double blood bag bucket (order no. 75007617)
- Centri-Lab® Type C round adapters
- Centri-Lab Type A rectangular adapters with the help of transition adapters that hold two Type A units per Thermo Scientific DoubleSpin™ bucket

LH-4000W rotor with double blood bag buckets



Safety, Reliability and Ease of Use in Blood Processing



Effortless access

The Soft-Touch Lid Lock ensures safe operation with an easy-open/close mechanism. The lid can be opened and closed effortlessly with just one hand.

Easy Set-Up, Programming and Monitoring

An intuitive operator keypad allows you to recall protocols at the touch of a button.

Whisper-Quiet Operation

A noise level under 58 dBA provides a more comfortable lab working environment.

Energy-Efficient Performance

Energy-saving design increases refrigeration efficiency and minimizes heat output, ensuring a more comfortable working environment.

Software for Sample Tracking

Heranet® software allows users to document centrifugation in blood banks. The software saves actual and set values and processes information which is entered at the centrifuge, using a barcode reader.

Highest Versatility

Swinging-bucket rotors for the Heraeus Cryofuge series centrifuges support a diverse array of sample processing applications without the need to modify existing protocols.

Our Heraeus Cryofuge 5500i centrifuge works with the Thermogenesis APX system from Thermogenesis company to automatically separate cord blood into a freezing bag, for the mononuclear cell product, an erythrocyte bag and an excess plasma bag.



SUPPLY
www.thermo.com/lab-pt

SPECIFICATIONS AND ORDERING INFORMATION

Thermo Scientific Heraeus Cryofuge 5500i Rotors and Adaptors

Rotor Specifications

Description	Capacity (place x mL)	Max. Tube Dim. D x L (mm)	Max. Speed (rpm)	Max. RCF (xg)	Cat. No.
LH-4000W Windshield Swinging Bucket Rotor	8 x 500	Blood Bags	4400	5480	75006476
LH-4000 Swinging Bucket Rotor	4 x 1000	One liter	4400	5350	75006475
Double Blood Bag Bucket (not applicable for 75006475)	8 x 500	Blood Bags	4400	5480	75006436
Round Bucket	4 x 1000	100 x 190	4400	5350	75006477
Biocontainment Lids (1 pair)					75006421
DoubleSpin™ Bucket	4 x 250	62 x 145	4400	4650	75006478
Biocontainment Lids (1 pair)					75006479
Microplate Adapter for DoubleSpin Bucket	4 x 8	Standard Microplates	4400	4650	75006486

Adaptors for LH-4000W and LH-4000 Rotors with Round Buckets

1.5 / 2 ml Disposable Microtubes	192 x 1.5/2 ml	11 x 50		75008132
7 ml Standard Tube	140 x 7 ml	13 x 177		75008133
5/7 ml Blood Collection Tube (Set of 4)	108 x 5/7 ml	13 x 177		11210822
10 ml Blood Collection Tube (Set of 4)	84 x 10 ml	16 x 177		11210821
10/15 ml Blood Collection Tube	68 x 10/15 ml	17 x 177		75008136
15 ml Standard Tube	76 x 15 ml	17 x 177		75008135
15 ml Disposable Conical Tube	48 x 15 ml	19 x 177		75008137
50 ml Round Bottomed Tube	28 x 50 ml	29 x 177		75006493
50 ml Disposable Conical Skirted Tube	20 x 50 ml	29 x 177		75006533
50 ml Disposable Conical Tube (Set of 4)	28 x 50 ml	29 x 177		11210826
225/175 ml Disposable Conical Bottle	4 x 225 ml	62 x 145		75008144 + 03095
250 ml Flat-Bottomed Bottle	4 x 250 ml	62 x 145		75008144
250 ml Disposable Conical Bottle	4 x 250 ml	62 x 177		75008147
500 ml Flat-Bottomed Bottle	4 x 500 ml	70 x 190		75008145
500 ml Disposable Conical Bottle	4 x 500 ml	96 x 190		75006438
1000 ml Flat Bottomed Bottle	4 x 1000 ml	98 x 190		no adapter needed

Adaptors for LH-4000W and LH-4000 Rotors with DoubleSpin buckets

1.5 / 2 ml Disposable Microtubes	448 x 1.5/2 ml*	11 x 50		75006452
7 ml Standard Tube	224 x 7 ml*	13 x 151		75006453
5/7 ml Blood Collection Tube	160 x 5/7 ml*	13 x 150		75006454
10/15 ml Blood Collection Tube	128 x 10/15 ml*	17 x 147		75006455
15 ml High Capacity Disposable Conical Tube	96 x 15 ml	19 x 150		75006468
15 ml Disposable Conical Tube	72 x 15 ml*	19 x 150		75006456
50 ml Round Bottomed Tube	32 x 50 ml*	29 x 140		75006491
50 ml Disposable Conical Tube	32 x 50 ml*	29 x 140		75006461
175/ 225 ml Disposable Conical Tub	e8 x 175/225 ml*	62 x 145		75006465 + 03095
250 ml Flat-Bottomed Bottle	8 x 250 ml*	62 x 145		75006465

*In order to achieve maximum specified capacity must order 8 adaptors (2/bucket)

Thermo Scientific Heraeus Cryofuge 5500i Rotors and Adaptors



Rotor Specifications					
Description	Capacity (place x mL)	Max. Tube Dim. D x L (mm)	Max. Speed (rpm)	Max. RCF (xg)	Cat. No.
Equipment for blood bag centrifugation in LH-4000W Rotor with Blood Bag Buckets					
Adapter for 500 ml quint, quad, triple with in-line filters and filter holder "XXL" (set of 6)	8 x 500 ml				76007677
Adapter for 450/500 ml quad, triple "XL" (set of 6)	8 x 500 ml				76007657
Adapter for 450 ml triple "L" (set of 6)	8 x 450 ml				76007647
Adapter for 450 ml double "M" (set of 6)	8 x 450 ml				76007667
Balancing weight system: taring range 300 – 1,550 g					75007668
Taring sticks: 4 of 6 g and 15 g					75007645
Taring plates: black rubber, 2 of 35 g and 2 of 65 g					75005759
Volume adjustment plates: blue rubber (set of 12)					75006681
Accessories for blood bag centrifugation in LH-4000/LH-4000W with round buckets					
Adapter for 400 – 500 ml quad blood bags "L/M" (set of 2)	4 x 500 ml				75006496
Adapter for small volume (250 ml) special blood applications (set of 2)	4 x 250 ml				75006485
Taring plates: black rubber, 2 of 35 g and 2 of 65 g					75005759
Taring sticks: 4 of 6 g and 15 g					75007645
Volume adjustment plates: blue rubber (set of 12)					75006681
BIOShield® Windshielded Rotor	4 x 250	62 x 125	5850	6960	75006435
Adaptors for BIOShield Rotor					
1.5/2 ml Disposable Microtubes	224 x 1.5/2 ml	11 x 50			75006452
7 ml Standard Tube	112 x 7 ml	13 x 117			75006453
5/7 ml Blood Collection Tube	80 x 5/7 ml	13 x 117			75006454
10/15 ml Blood Collection Tube	64 x 10/15 ml	17 x 117			75006455
15 ml Disposable Conical Tube	36 x 15 ml	19 x 119			75006456
50 ml Round Bottomed Tube	16 x 50 ml	29 x 119			75006491
50 ml Disposable Conical Tube	16 x 50 ml	29 x 119			75006461
175/225 ml Disposable Conical Tube	4 x 175/225 ml	62 x 125			75006465 + 03095
250 ml Flat-Bottomed Bottle	4 x 250 ml	62 x 125			75006465
HIGHplate® Windshield Rotor	2 x 5 Standard Microplates		5650	5890	75006435
Standard Microplates	10 x Standard Plates up to 20 mm tall				
Deepwell Microplates	2 x Deepwell Plates up to 83 mm tall				
Diagnostik™ Rack Rotor			3500	2550	75006480
Trays for Diagnostik Rotor					
Hitachi	100	118 x 20			75006416
Olympus	120	117 x 20			75006417
Symex	120	200 x 25			75006418
Dade-Behring	40				75006419



Rotor Specifications

Description	Capacity (place x mL)	Max. Tube Dim. D x L (mm)	Max. Speed (rpm)	Max. RCF (xg)	Cat. No.
 Highconic® fixed-angle Rotor	6 x 50	29 x 119	8500	10,020	75003057
Adaptors for Highconic rotor					
1.5 ml Disposable Microtubes	24 x 1.5 ml	11 x 58			75002905
3.5 ml Round Bottomed Tube	24 x 3.5 ml	11 x 103			75003091
7 ml Round Bottomed Tube	12 x 7 ml	13 x 115			75003092
12 ml Round Bottomed Tube/5 ml Blood Collection Tube	12 x 12/5 ml	16 x 96			75003093
15 ml Disposable Conical Tube	6 x 15 ml	17 x 119			75003095
16 ml Round Bottomed Tube	6 x 16 ml	18 x 124			76002906
38 ml Round Bottomed Tube	6 x 38 ml	25 x 112			75003094
50 ml Round Bottomed Tube	6 x 50 ml	29 x 118			75003102
50 ml Disposable Conical Tube	6 x 50 ml	29 x 119			75003103
LAC-250 Fixed-angle rotor	6 x 250	62 x 145	10,000	15,320	75006483
Adaptors for LAC-250 rotor					
3 ml Glass Tube	30 x 3 ml	10 x 75			00370 + 00456
5 ml Glass Tube	72 x 5 ml	12 x 75			00388
10 ml Round Bottomed Tube	42 x 10 ml	16 x 100			00389
15 ml Round Bottomed Tube	30 x 15 ml	18 x 120			00456
30 ml Round Bottomed Tube	18 x 30 ml	24 x 105			00449
50 ml Disposable Conical Tube	6 x 50 ml	29 x 119			03072
175/225 ml Disposable Conical Bottle	6 x 175/225 ml	62 x 130			03095
250 ml Flat-Bottomed Bottle	6 x 250 ml	62 x 130			N/A





Thermo Scientific Heraeus Cryofuge 5500i

Specifications		
Description	Speed range rpm 300 – 10,000, adjustable in 10 rpm increments	
Max. RCF	15,320 g	
Max. capacity	8 x 500 ml blood bags or 4 x 1,000 ml or 4 x 8 standard microplates	
Drive	brushless induction drive	
Accel./decel. profiles	9/9	
Program memory	9 + 1 pretemp program	
Run time	0 – 99 min, plus hold	
Temperature range	-9 to +40 °C , CFC-free refrigerant	
Control	microprocessor controlled via EASYcontrol™ II	
Functions	RCF-preselection, Quick-run, automatic rotor recognition, imbalance detection and Soft-Touch lid-lock	
Imbalance	electronic, speed and rotor dependent	
Dimensions (H x W x D)	mm	960 x 680 x 720
	inch	37.8 x 26.8 x 28.3
Weight	kg (lbs)	286 (631)
Heraeus Multifuge 4KR		
Dimensions (H x W x D)	mm	700 x 680 x 720
	inch	37.8 x 26.8 x 28.3
Heraeus Multifuge 4KR weight	kg (lbs)	236 (519)
Test standards		
Designed and tested in accordance with EN 61 010-1, EN 61 010-2-020, EN 50 081-1, EN 50 082-1		

Ordering Information		
Model	Description	Order no.
Heraeus Cryofuge 5500i	Refrigerated, 230 V, 50 Hz, 12.9 A	75004471
Heraeus Cryofuge 5500i (CSA)	Refrigerated, 200 – 230 V, 60 Hz, 16 A / 200 V, 50 Hz, 16 A	75004479
Heraeus Cryofuge 5500i (CSA)	Refrigerated, 120 V, 60 Hz, 20A	75004475
Heraeus Multifuge 4KR ¹⁾	Refrigerated, 230 V, 50 Hz, 12.9 A	75004461
GMP Version		
Heraeus Cryofuge 5500i (water cooled)	230 V/50 Hz /single phase	75004478

Heraeus Cryofuge 5500i GMP with complete encapsulated housing and special smooth surface meets requirements for "clean rooms" according to EC-GMP Guide 1 Annex and FDA Aseptic Guide

¹⁾ Not available in the United States

If floor space is limited, consider our Heraeus Multifuge 4KR. Offering the same specifications as our Cryofuge 5500i, but sized to fit under your lab workstation.



© 2008 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

www.thermo.com/cryofuge

North America: USA/Canada +1 866 984 3766

Europe: Austria +43 1 801 40 0, Belgium +32 2 482 30 30, France +33 2 2803 2000, Germany national toll free 08001-536 376, Germany international +49 6184 90 6940, Italy +39 02 02 95059 434-254, Netherlands +31 76 571 4440, Nordic countries +358 9 329 100, Russia/CIS +7 (812) 703 42 15, Spain/Portugal +34 93 223 09 18, Switzerland +41 44 454 12 12, UK/Ireland +44 870 609 9203

Asia: China +86 21 6865 4588 or +86 10 8419 3588, India +91 22 6716 2200, Japan +81 45 453 9220, Other Asian countries +852 2885 4613

Countries not listed: +49 6184 90 6940 or +33 2 2803 2000